



# FP 250 DIGITAL

## PROFESSIONAL COMPRESSION THERAPY



### DEVICE AND ITS GARMENTS

Digital base unit	250UBD
Left leg sleeve/7 sectors	250G7L
Right leg sleeve/7 sectors	250G7R
Left foot sleeve	000PNVL
Right foot sleeve	000PNVR
Arm Sleeve (8 sectors)	250B8
Arm sleeve with separate Hand (8 sectors)	250B6M2
Abdominal strip (2 sectors)	000FA2
Professional Duplication Set	100SD8

## USE AND INDICATION

The Fisiopress FP 250 DIGITAL device is an outpatient system highly recommended for venous and lymphatic pathologies, to drain the liquid component of fluids, deflate the limb and restore normal pressure balance. The operative cycle, rapid and efficient, together with the practical and ergonomic design of the garments (easily manageable by the operator via the simple professional duplicating set) allows you to exercise effective compression and reduce edema in a short time. The system, despite being professional is small in size and can be easily placed in all medical surgeries/ or clinics and does not require any maintenance. It can be easily moved together with the group of garments without removing any connection, thus facilitating the management of the medical bed when it is to be used for different therapies.


In professional use they are available simple guidelines for an appropriate use of the instrument, which, related to an appropriate scientific literature, allow to obtain the best results of this method. It is always advisable to include the compression therapy within a protocol integrated therapeutic.

This device can also drive an arm sleeve with a separate hand for the treatment of secondary lymphedema particularly present in the patient's hand

The device has a medical CE mark



ARM WITH SEPARATE HAND

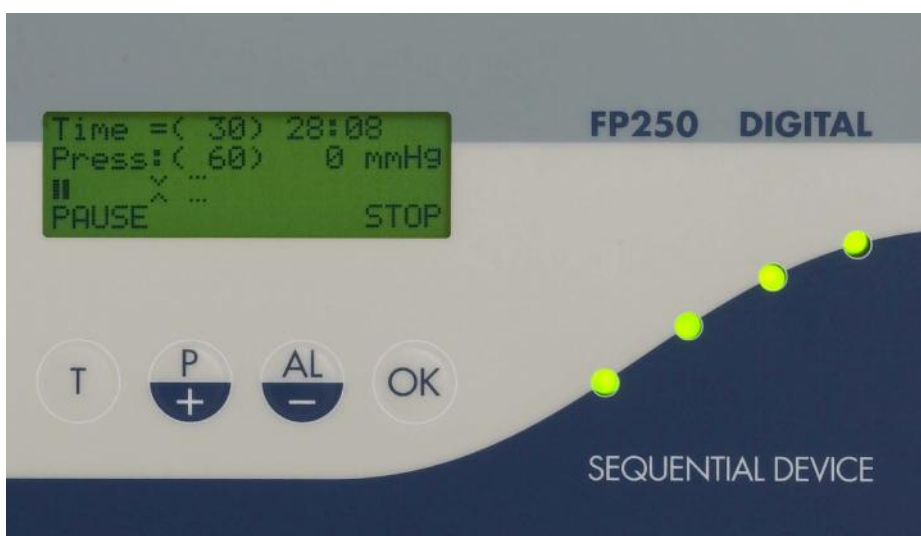

DIGITAL DISPLAY



INTERNAL PNEUMATIC BAGS PARTIALLY OVERLAPPED

## TECHNICAL FEATURES

- ) Sequential compression therapy with distal-proximal pressure gradient distributed on 8 outlets, capable to simultaneously treat one or two limbs.
- ) Pressure range: 0 - 160 mm of Hg. (on an applied garment).
- ) Operating cycle: 30 seconds.
- ) Large digital display for viewing pressure and time values programmed treatment and remaining time.
- ) Pre-set alarm thresholds.
- ) Graph of pressures and programmable functions.
- ) High pressure alarms (audible and automatic shutdown).
- ) Alarms on low pressures (acoustic, due to disconnection of connection..)
- ) Electronic adjustment of the treatment time, both automatically and manually.
- ) DASC system (Distribution to communicating sectors): this system, unique to Fisiopress, allows have an always homogeneous pressure in all the chambers as they fill up with air, thus allowing the terminals to adapt to the geometry of the treated limb and compress it evenly without danger of stagnation or reverse flows.
- ) Automatic shutdown at the end of therapy during deflation.
- ) Garments available for lower limb: leg and foot (8 sectors), separated for a better fit and adaptation to the length of the treated limb. Four possibilities of size adjustment circulate through special zip fasteners.
- ) Possibility of connecting a couple of leg sleeves and the abdominal strip at the same time.
- ) 2 Sectors Abdominal strip (2 sectors)
- ) Garments available for upper limb: Arm Sleeve (8 sectors). Three adjustment possibilities of circular size via special zip fasteners.
- ) All garments have parallel and partially overlapped internal pneumatic sectors.
- ) Possibility of individual replacement of each internal pneumatic bag.
- ) Entirely hand-sewn garments.
- ) Medical CE mark