



## FP 4 SSP (Slow Sequential Press)

SEQUENTIAL COMPRESSION THERAPY  
BEAUTY & WELLNESS





## DEVICE AND GARMENTS

Digital Base Unit	004SSPDIGUB
Leg sleeve (6 sectors) Right	004SSPDIG6R
Le sleeve (6 sectors) Left	004SSPDIG6L
Arm sleeve (8 settori)	004SSPDIGB8
Abdominal strip (2 settori)	004SSPDIGFA2
Professional Duplication Set	004SSPDIGSD8


## USE AND APPLICATIONS

The top-of-the-range professional device, **Fisiopress FP 4 SSP** with 8 output sectors, in the Digital version (in the photo) is the best sequential technology allows today on this number of sectors. This system is indicated:

- Treatment of water retention.
- The treatment of cellulite.
- Treatment of venous insufficiency and all its consequences (swollen and painful legs, ulcers, oedema, varicose veins, etc...).
- Lymphostasis therapy (lymphoedema, edema, etc...)
- Prevention of such problems and their worsening.
- Possibility to combine with other aesthetic treatments or medical therapies used to prevent or reduce the blemishes and pathologies described above.



The philosophy of the Fisiopress FP 4 SSP system follows a setting defined as "natural pressotherapy" in how much the massage it performs pushes the fluids that stagnate along their "natural" direction, an action entirely congruous with the physiological one of our organism. The concept of natural pressure therapy is therefore based on an action performed by the most similar device possible to that exercised by a good operator who sets himself the goal of discharging the liquids that stagnate inside the limb to be treated, putting the physical mechanism back in motion, it is no longer well functioning by pushing the fluids in their physiological direction. The 8 sectors in which leg sleeves/abdominal strip and arm sleeve (end pieces) are divided allow a continuity and homogeneity of the thrust that no other pressure therapy device is capable of exercise. Not only is the action generated highly effective, but the sensation felt by the person subjected to the treatment is very pleasant and extremely relaxing.

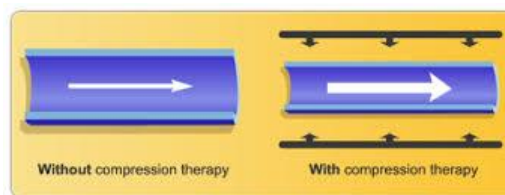

## RESULTS

The results obtainable using a truly sequential pressure therapy system with natural action are:

- drainage of excess liquids
- restoration of water exchanges
- decrease in limb volume
- ) decrease in pain
- prevention of cellulite
- feeling of lightness in the limb
- improvement of venous circulation
- improvement of lymphatic circulation



## TECHNICAL FEATURES

- Sequential compression therapy with distal-proximal pressure gradient distributed over 8 outlets, capable of treating one or two limbs simultaneously and abdominal fascia.
- Pressure range: 0 - 100 mm of Hg
- Operating cycle: 120 seconds.
- Large backlit digital display for viewing pressure values, programmed treatment time and residual time, and pre-set alarm thresholds.
- Graph of pressures and programmable functions.
- High pressure alarms (audible and automatic shutdown.
- Alarms on low pressures (acoustic, due to disconnection of connection..)
- Electronic adjustment of treatment time, both automatically and manually.
- D.A.S.C. system (Distribution to communicating sectors): this system, unique to Fisiopress, allows for an always homogeneous pressure in all the chambers as they fill up with air, thus allowing the terminals to adapt to the geometry of the treated limb and compress it uniformly without danger of stagnation or reverse flows.
- Automatic shutdown at the end of therapy during deflation.
- Terminals available for the lower limb: leg and foot (6 sectors),
- Possibility of connecting a pair of leggings and the abdominal band in 2 sectors at the same time.
- Garments available for upper limb: arm sleeve (8 sectors). Three possibilities of adjusting the circular size using special zip fasteners.
- All garments have parallel and partially overlapping internal sectors.
- Possibility of individual replacement of each internal pneumatic bag.
- Entirely hand-sewn garments.
- Aesthetic certification

